



a. MCDP (10 - *N* - methylcarbamoyl - 3, 7 - dimethylamino - 10 - *H* - phenothiadine)

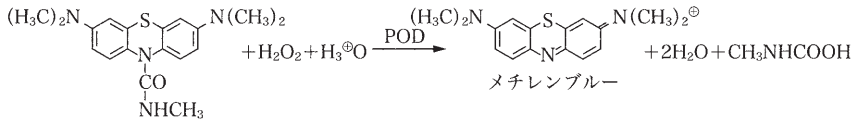

b. BCMA [Bis (3 - bis (4 - chlorophenyl) methyl - 4 - dimethyl - aminophenyl) amine]

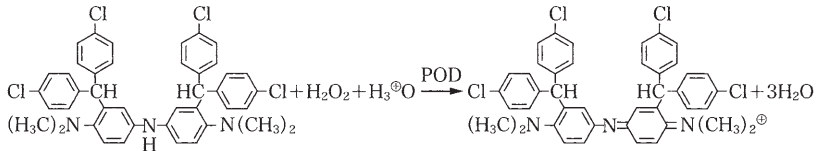

	生成色素		感度比 (Phenol = 1)	検査試薬項目(応用例)
	$\lambda \text{ max}$	$\epsilon \times (10^4)$		
MCDP	666	9.6	18.2	Dupan-2, HTLV-Iなど
BCMA	755	6.3	11.9	NEFA, UA, NAGなど

図 6-1 色素酸化反応法におけるロイコ型発色試薬の例